

Molift Duo

A free standing and flexible gantry system

oetac[®]

Approved installer





Molift Duo

With a Molift free standing unit, we offer a system that is easy to install without additional fixtures on ceiling or walls. It is ideally suited to temporary lift applications in hospitals, nursing homes or private homes.

The easy solution for all bed and room sizes

Molift Duo offers a system that is easy to install without being fixed to the ceiling or walls. It is ideally suited to temporary hoist applications in hospitals, nursing homes or private homes. In combination with Molift Air or Molift Nomad hoist and a sling, Molift Duo solves most overhead lifting situations, whatever the environment. With a lifting capacity of 300 kg, it is ideal for lifting and transferring a person from or to; bed, floor, chair, wheelchair, toilet etc.

Lightweight construction

Molift Duo has a lightweight aluminium construction and a functional design which makes the gantry system easy to install, dismantle and move elsewhere.

This 2-post gantry system has a rail length between 2000 mm and 3500 mm and has 4 step adjustable heights up to 2550 mm, making this free standing hoist unit suitable for all bed and room sizes.



| Product name | | Item no. | Weight | | |
|--|------------------------------|--|---|----------------------|----|
| Molift Duo 2000 Molift Duo 2500 Molift Duo 3000 Molift Duo 3500 | | 1600013 1600015 1600017 1600019 | Maximum user weight (SWL) Total weight Material Anodized aluminium, Plastic | 300 kg 38 - 44 kg | |
| Weight of unit Post (complete) Post (without foot) Boom | 14,8 kg 11 kg 4,2 kg/m | | Alk. pH 300 kg | 661 lbs | CE |
| Foot | 3,8 kg | | Clean / disinfect the product using a standard non-abrasive cleaning agent with a pH between 5 and 9, or with a 70% disinfectant solution. The product can be safely cleaned in a cabinet washer-disinfector at 85 °C for 3 minutes. | | |



Features



Easy to install

Due to the lightweight system, it is easy to install the gantry and move it to a new location.



Supreme stability

The gantry system has supreme stability and a functional design.




Adjustable height


Height positions can be set at four fixed positions from 225 cm up to 255 cm. The height is adjustable with the quick release pin.



Charging of lift motor

Molift hoist motors can easily be charged thorough the hand control or optionally through in-rail charging.

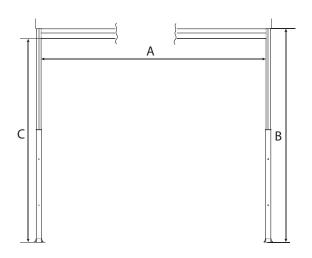

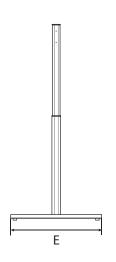

Measurement

 A: Rail lengths
 2000, 2500, 3000, 3500 mm

 B: Height
 2250 - 2550 mm

 C: Internal height
 2400 mm

 E: Depth
 1200 mm

 Max internal space
 3500 x 2400 mm





Etac is a world-leading developer and provider of ergonomic assistive devices and patient handling equipment. Our heart lies in the solutions that optimise quality of life for the individual, their family and caregivers.



